

PHOLA PRO POSITIONING WHEELCHAIR

PRODUCT CODE: 12376



The Phola Pro Tilt-in-Space Wheelchair is engineered for advanced positioning, comfort, and durability with folding capabilities. Its lightweight aluminium folding frame combined with a tension backrest provides supportive and comfortable seating. The chair features 25° seat tilt and up to 160° reclining, allowing effective pressure relief and positioning throughout the day. Flip-back detachable armrests, elevating leg rests, and adjustable push handles allow easy transfers and caregiver convenience. A foldable handle brake, adjustable headrest, and optional stabilization bar enhance safety, control, and user support.



TAB



Headrest



Elevating Leg Rests

Product Specifications




Recline:
160 degrees



Tilt:
25 degrees



Seat Width:
16", 18", 20"

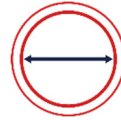

Rear Wheel:
24" Spoke
Wheel



Front Wheel:
7" Castor




Seat Height:
51cm



Length:
114cm



Total Width:
Seat Width +25cm

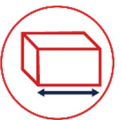

Packaging DIM:
73 x 22.5 x 86cm



Max Weight Capacity:
100 / 125kg



Net Weight:
33kg



Gross Weight:
45kg

Medop

www.medop.co.za | [@medopcc](https://www.instagram.com/medopcc)



Johannesburg:
011 827 5893
enquiries@medop.co.za
Medopark, 1 Manchester
Road, Wadeville

Cape Town:
021 534 8628
admin@medop.co.za
Viking Business Park, Corner
Park & Viking Road, Epping