


# ZAMBA KNEE WALKER

PRODUCT CODE: HCT9150



The Zamba Knee Walker is built for comfort and convenience during recovery from foot, ankle, or lower leg injuries. Its durable steel frame and automotive-style tie-rod steering provide stability and control. The contoured knee platform fits either leg, while the 10" front and 8" rear TPU wheels offer smooth, all-terrain mobility. For easy transport, the handlebar and rear wheel detach quickly.

With tool-free adjustments, the handlebar height ranges from 31.1" to 38.1" and the knee pad from 14.6" to 18.7", allowing a personalized fit. The cushioned 12.4" x 6.3" knee pad offers day-long support and comfort. Compact, adjustable, and easy to manoeuvre, the Zamba Knee Walker helps you stay mobile and independent during healing.

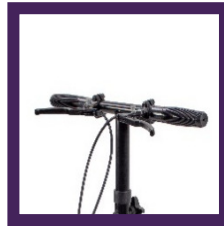

Automotive-style tie-rod steering



Adjustable knee pad  
31x17cm



10" Front Castor


## Product Specifications



Total Height:  
79cm – 96cm



Rear Wheel:  
8"

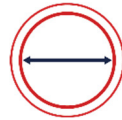

Front Wheel:  
10"



Knee Pad Height:  
41cm – 49cm



Packaging DIM:  
46.5 x 42.5 x  
22cm



Max Weight Capacity:  
105kg



Net Weight:  
8.9kg



Gross Weight:  
10kg

[www.medop.co.za](http://www.medop.co.za) | [@medopcc](https://www.instagram.com/medopcc)



**Johannesburg:**  
011 827 5893  
[enquiries@medop.co.za](mailto:enquiries@medop.co.za)  
Medopark, 1 Manchester Road, Wadeville

**Cape Town:**  
021 534 8628  
[admin@medop.co.za](mailto:admin@medop.co.za)  
Viking Business Park, Corner Park & Viking Road, Epping